


Determine the asymmetry factor in HPLC at 10 percent peak height - How To

Measure the left half of the peak width at 10% peak height (A) and then measure the right half of the peak width at 10% peak height (B) (see figure below). The ratio of B over A is the asymmetry factor. A value of 1.0 is perfect symmetry and 0.9 – 1.2 is acceptable.

Value less than 1 is for a fronting peak and greater than 1 is for a tailing peak. Values greater than 2 are unacceptable.

ASYMMETRY FACTOR

$$A_f = \frac{B_{(10\% h)}}{A_{(10\% h)}}$$

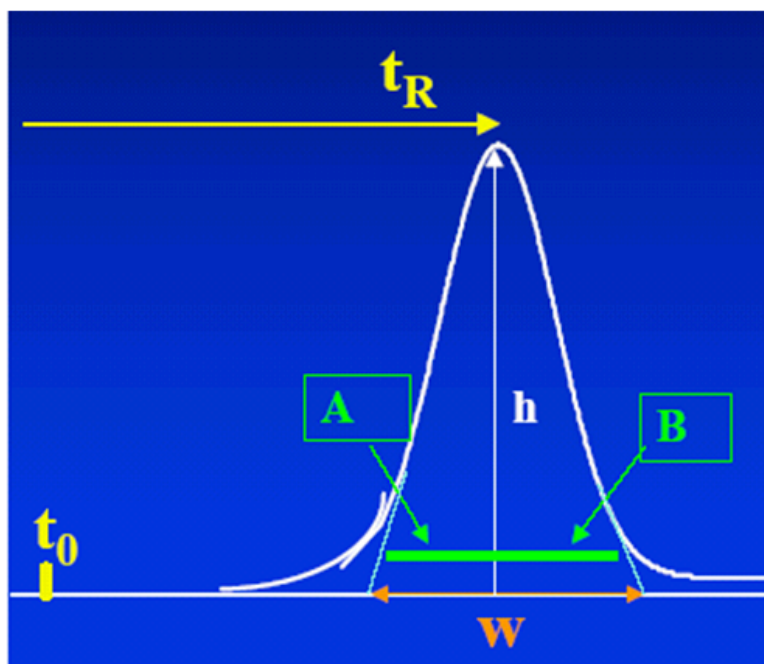

COGENT 
HPLC Columns™