


## Polyphenols Gallic Acid, Catechin, and Epicatechin in Grape Seed Extract Analyzed with LCMS - AppNote

### Separation of Gallic Acid, Catechin, and Epicatechin

This Method was applied to the determination of Phenolics in a grape seed extract. From *Figure B*, it can be determined that the Aqueous extract contains a high amount of Gallic Acid, Catechin and Epicatechin.

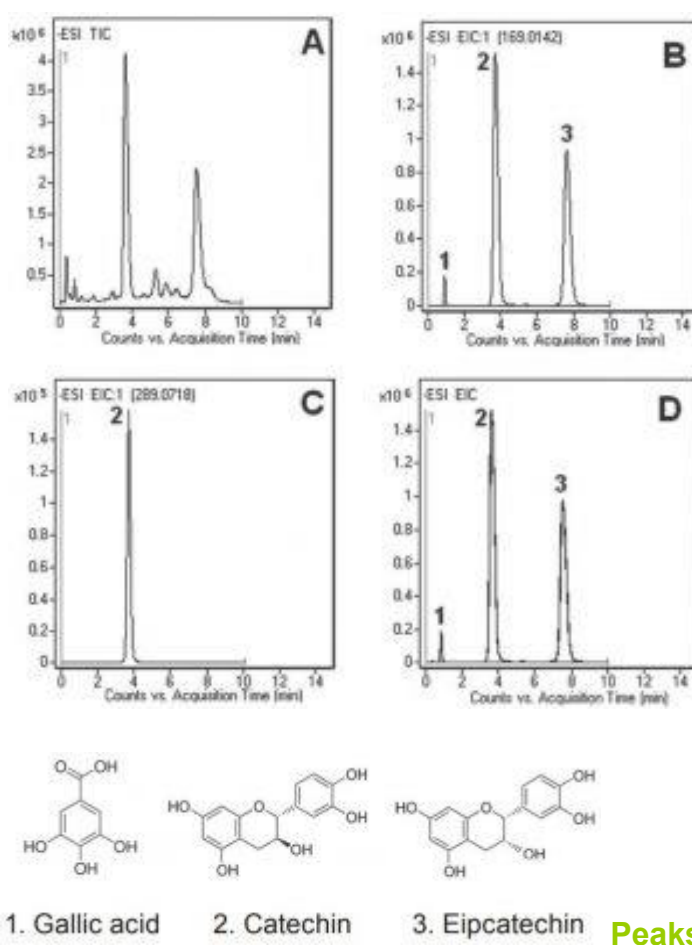

**Peaks:**

1. Gallic Acid 169.0142 m/z (M-H)<sup>-</sup>, RT=0.862 min
2. Catechin 289.0718 m/z (M-H)<sup>-</sup>, RT=3.739 min
3. Epicatechin 289.0718 m/z (M-H)<sup>-</sup>, RT=7.626 min

### Method Conditions

**Column:** Cogent Bidentate C18™, 4μm, 100Å

**Catalog No.:** [40018-05P-2](#)

**Dimensions:** 2.1 x 50mm

### Mobile Phase:

A: DI Water / 0.1% Formic Acid

B: Methanol / 0.1% Formic Acid

**Gradient:**

Time (minutes)	%B
0	10
5	15
7	15
8	10

**Injection vol.:** 1µL

**Flow rate:** 0.4mL / minute

**Detection:** ESI – NEG - Agilent 6210 MSD TOF Mass Spectrometer

**Sample Preparation:** Grape seed extract tablet was crushed and extracted with DI Water at 40°C. Before injection, sample solution was filtered using a 0.45µm Nylon Syringe Filter (MICROSOLV Tech Corp.).

**t<sub>o</sub>:** 0.4 minutes

---

**Note:** Grape seed extracts have a high content of phenolic compounds, such as Gallic Acid, Catechin, Epicatechin and Procyanidins. They have been the matter of intense investigations with respect to their potentially beneficial effects on human health. Phenolic compounds are secondary plant metabolites and they play an important role in plant growth, reproduction, and protection against pathogens and predators. This data was presented at the American Society of Pharmacognosy Annual Meeting and Exhibition, July 30 – August 3, 2011, San Diego, California.

---



**Attachment No 138 Polyphenols in Grape Seed Extract Analyzed with LCMS pdf 0.3 Mb**

[Download File](#)

Printed from the Chrom Resource Center  
Copyright 2025, All Rights Apply  
**MicroSolv Technology Corporation**  
9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900  
Fax: (910) 769-9435  
Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)  
Website: [www.mtc-usa.com](http://www.mtc-usa.com)