
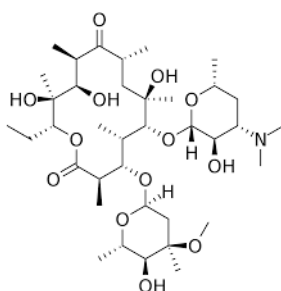


Erythromycin Analyzed with HPLC ELSD - AppNote

Erythromycin, a macrolide antibiotic

As this compound lacks chromophores, it typically requires derivatization for use in UV detection. In this method, we retain excellent peak shape without the need for these pre-column derivatization steps saving time and resources.

RSD values (less than 0.5%) demonstrate the consistent and reliable retention, as shown in the 10 injections overlay below.

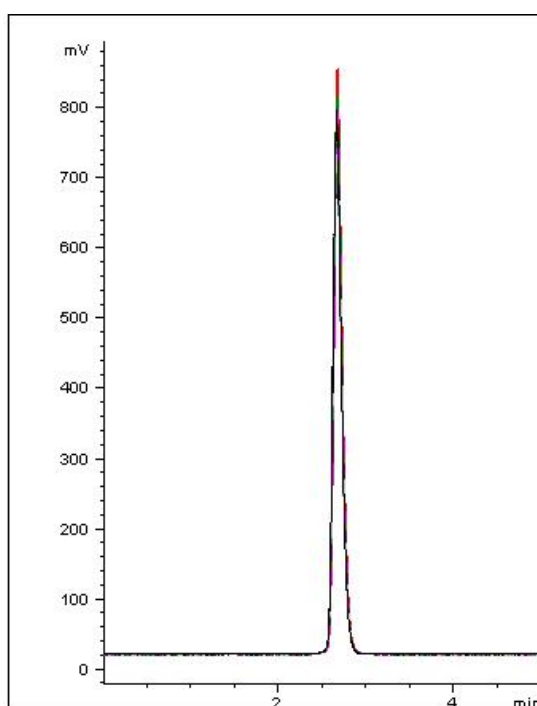

Peak:

Erythromycin

Method Conditions:

Column: Cogent Diamond Hydride™, 4µm, 100Å.

Catalog No.: [70000-10P](#)

Dimensions: 4.6mm x 100mm

Mobile Phase: 50% acetonitrile / 50% DI water / 0.1% formic acid

Injection vol.: 1 µL

Flow rate: 1.0mL / minute.

Detection: ELSD, gain: 9, temperature: 50°C, nitrogen: 3.5 bar.

Sample Preparation: 2.0mg / mL Erythromycin HCL in DI water.

Notes: Erythromycin is a macrolide antibiotic used to treat bacterial infections and is also used to prevent recurrent rheumatic fever.



Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com