
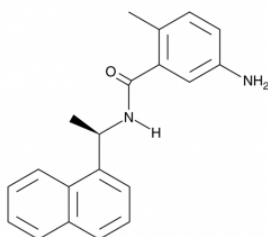


## Covid Drug Molecule GRL0617 Assay with HPLC - AppNote

### GRL-0617, a benzamide anti-Covid molecule and potent non-covalent inhibitor

This simple gradient analyzing GRL-0617 is easily performed, and the 5 injections overlay below demonstrates run to run consistency and great peak symmetry. The peak efficiency is very good, peak tailing is minimal, and repeatability is easily met with RSD values being less than 1.0%.

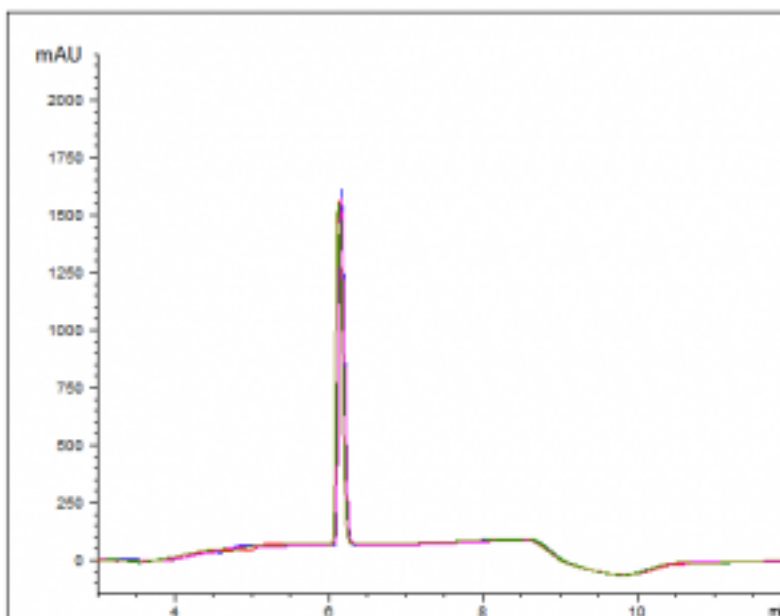

#### Peak:

GRL0617

#### Met hod Conditions:

**Column:** Cogent RP Phenyl Hexyl™, 5um, 100Å

**Catalog No.:** [68539-15P](#)

**Dimensions:** 4.6mm x 150mm

#### Mobile phase:

A: DI water and 0.1% formic acid

B: Acetonitrile and 0.1% formic acid

## Gradient:

time (minutes)	%B
0	25
1	25
5	85
6	85
8	25
10	25

**Flow rate:** 1.0mL / minute

**Injection volume:** 1 µL

**Detection:** UV @ 224nm

**Sample Preparation:** 1.0 mg/mL of GRL-0617 in methanol.

**%RSD:** <0.5%

**t<sub>0</sub>:** 2.0 minutes

**K':** 2.075

---

**Note 1:** GRL0617 is an inhibitor of SARS-CoV and SARS-CoV-2 papain-like proteases (PL-pro) . It is selective for the viral proteases over human proteases. GRL0617 blocks both PL-pro-mediated actions towards host biological processes and virus replication.

**Note 2:** Capacity is determined using the following equation:  $k = (t_R - t_0)/t_0$

- $t_R$  = Retention Time of an Analyte Peak
- $t_0$  = Retention Time of non-Retained Peak



Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)