
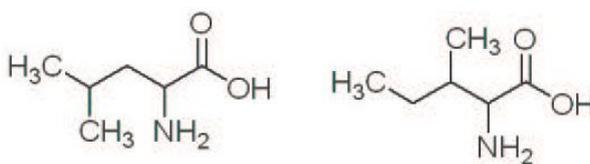


## Leucine and iso Leucine Separated with LCMS - AppNote

### Leucine and Isoleucine, Isobaric Compounds Separated for LCMS Detection

This Application Note presents a Method for the Analysis of two Branch Chain Amino Acids (BCAA) with nearly identical  $m/z$  values.

Addition of weak Acetic Acid to the Mobile Phase allows the small difference in the pKas of the two Isobaric Compounds to achieve Resolution. A few microliters of Serum suffices as the Sample to estimate the Leucine and Isoleucine present. The two compounds are Separated and can be easily Quantified.

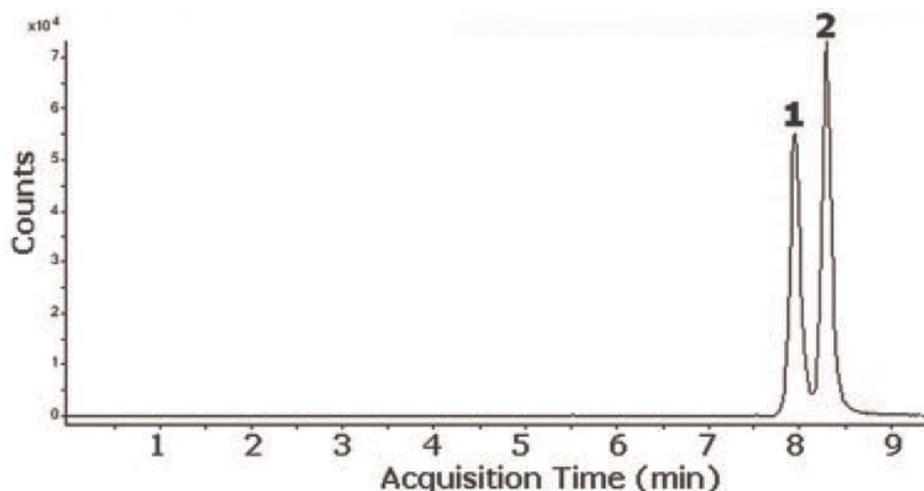

#### Peaks:

1. D,L - Leucine 132.1025  $m/z$  (M+H)<sup>+</sup>,
2. D,L - Isoleucine 132.1025  $m/z$  (M+H)<sup>+</sup>

#### Method Conditions

**Column:** Cogent Diamond Hydride™, 4  $\mu$ m, 100 Å

**Catalog No.:** [70000-15P-2](#)

**Dimensions:** 2.1 x 150 mm

#### Mobile Phase:

A: DI Water / 0.1% Acetic Acid

B: Acetonitrile / 0.1% Acetic Acid

**Gradient:**

Time (minutes)	%B
0-2	85
2-3	80
3-5	80
5-6	75
6-7	75
7-8	70
8-9	70
9-10	50
10-11	50
11-11.1	20
11.1-14	20

**Post Time:** 5 minutes**Flow rate:** 0.4 mL / minute**Detection:** ESI – pos - Agilent 6210 MSD TOF Mass Spectrometer**Sample Preparation:** 100 ng / mL of each Compound prepared in 50:50 Solvent A / Solvent B and diluted 1:10.

**Notes:** D,L - Leucine and D,L - Isoleucine are considered essential Amino Acids because human beings cannot survive unless these Amino Acids are present in the diet. The Method is especially applicable to the frequent monitoring of Serum while treating patients with Maple-Syrup Disease, and may also be used for rapid diagnosis of the disease in suspected infants.

[1] Chromatogram is adapted from: J.J.Pesek, M.T. Matyska, S. M. Fisher, T. R. Sana, *Journal of Chromatography A*, 1204 (2008) p50.



**Attachment No 84 Leucine & iso Leucine Separated with LCMS pdf** 0.1 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)Website: [www.mtc-usa.com](http://www.mtc-usa.com)