



## Fexofenadine Analyzed with HPLC - AppNote

### Simple Isocratic Assay Method without Amine Additives

Click [HERE](#) for Column Ordering Information.

Fexofenadine, antihistamine, has a number of aromatic moieties and the tertiary amine in its chemical structure that can cause Peak tailing issues when using most L11 Columns for Analysis. The USP Assay Monograph for Fexofenadine Tablets uses Triethylamine to reduce tailing.

With this Method, Trifluoroacetic Acid is used instead and good Peak shapes are observed, as the five Chromatograms that are overlaid in the Figure show.

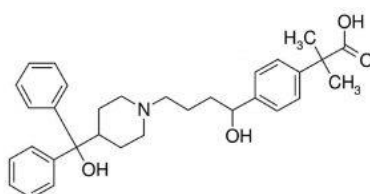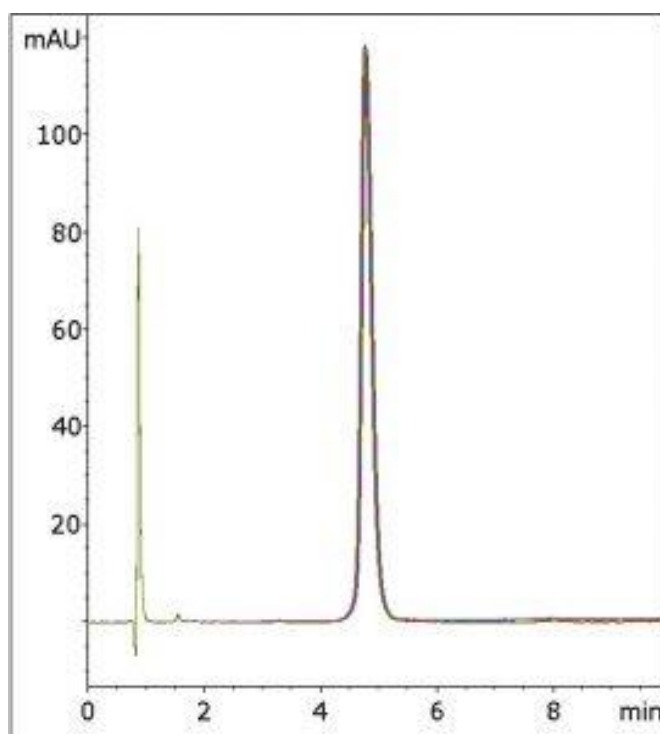

**Peak:**

Fexofenadine

### Method Conditions

**Column:** Cogent Phenyl Hydride™, 4μm, 100Å

**Catalog No.:** 69020-7.5P

**Dimensions:** 4.6 x 75mm

**Mobile Phase:** 62% DI Water / 38% Acetonitrile / 0.1% Trifluoroacetic Acid (TFA)

**Temperature:** 35°C

**Flow rate:** 1.0mL / minute

**Detection:** UV @ 220nm

**Injection vol.:** 1µL

**Sample Preparation:**

*Stock Solution:* 180mg strength Allegra® (Fexofenadine) tablet was ground and added to a 100mL volumetric flask. It was diluted with the Mobile Phase and vortexed 5 minutes. A portion was filtered with a 0.45µm Nylon Syringe Filter (MICROSOLV Tech Corp.).

*Working Solution:* A 100µL aliquot of the *Stock Solution* was diluted with 900µL of the Mobile Phase.

**t<sub>o</sub>:** 0.87 minutes

---

**Note:** Fexofenadine is marketed under the trade name Allegra, but generic versions are also available. In January 2011, the Food and Drug Administration approved over-the-counter sales of the drug without a prescription. It is a widely selling antihistamine used for treatment of hay fever and other allergies. Because it does not cross the blood-brain barrier, it causes less drowsiness than first generation antihistamines.

---



**Attachment No 139 Fexofenadine Analyzed with HPLC pdf** 0.3 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)