
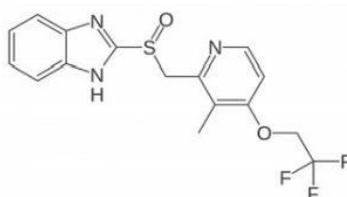


Lansoprazole Capsule Analyzed with HPLC - AppNote

Separation from matrix peaks

Lansoprazole is separated from several Matrix Components in this simple assay method. The Mobile Phase is LC-MS compatible so the method could be applied to more complex samples such as plasma. Five runs are shown in the figure to illustrate the repeatability of the data.

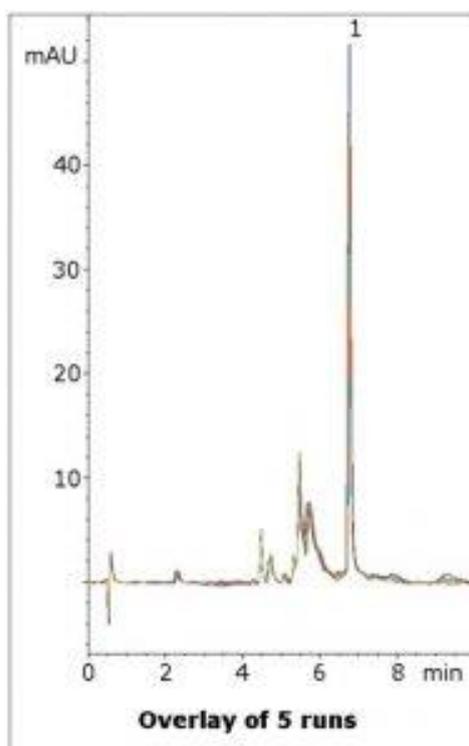

Peak:
Lansoprazole

Method Conditions

Column: Cogent Bidentate C18™, 2.2 μm, 120 Å

Catalog No.: [40218-05P-2](#)

Dimensions: 2.1 x 50 mm

Mobile Phase:

A: DI Water / 0.1% Formic Acid (v/v)

B: Acetonitrile / 0.1% Formic Acid (v/v)

Gradient:

Time (minutes)	%B
0	10
1	10
6	70
7	10

Injection vol.: 0.2 µL

Flow rate: 0.4mL / minute

Detection: UV @ 285 nm

Sample Preparation: 15 mg strength Lansoprazole capsule contents were ground and added to a 25 mL volumetric flask. A portion of 50/50 Solvent A / Solvent B was added and the flask was sonicated 10 minutes. Then it was diluted to mark and mixed. A portion was filtered through a 0.45 µm Nylon Syringe Filter (MICROSOLV Tech Corp.).

t_o: 0.5 minutes

Note: Lansoprazole is a proton-pump inhibitor used for acid reducing effects. It is available over-the-counter by Novartis under the trade name Prevacid, as well as generic versions.



Attachment No 273 Lansoprazole Capsule Analyzed by HPLC pdf 0.3 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com