


HPLC Pre-Column Filters

Low Cost, Effective, Safe and Disposable

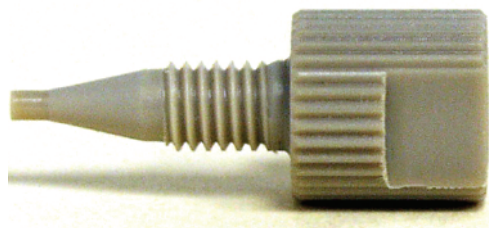

Protect Columns and
Minimize Troubleshooting Time

These HPLC “column savers” are designed to give you the very best filter capacity for un-inhibited flow while providing the best protection for your HPLC Columns. These easy to use filters prevent un-dissolved particles from getting into your HPLC column’s frits or into the column.

Hand Tighten – Knurled for easy and comfortable hand tightening and has a “universal” fitting so it will not add any system or dead volume.

Wrench Tighten – A “flat” area is provided for that “little extra” that is sometimes needed to make sure the connection is perfect.

Peek Construction – Compatible with most solvents for HPLC; these filters are durable.

Stainless Steel Filter – Made from woven stainless steel cloth these filter units provide high mechanical strength and are made from the same material as your HPLC columns.

Disposable – When your back pressure increases, the first step to take will be to replace these low cost filters with a fresh unit. If the back pressure returns to normal, you have solved the problem.



Easy to Use and install, simply place the pre-column filters inline as shown above, between your male nut and your column or guard column if you are using one.

Cat. No.	Description
49650-005	HPLC Pre-Column Filters with 0.5um SS Filter Elements. Universal Connection, Hand or Wrench Tightened. Disposable. 1 Each.
49650-020	HPLC Pre-Column Filters with 2.0um SS Filter Elements. Universal Connection, Hand or Wrench Tightened. Disposable. 1 Each.

L-64



MTC • MicroSolv Technology Corporation

World Headquarters

1 Industrial Way West, Bldg E, Eatontown, NJ 07724 USA

Ph. 1-732-578-1777 Fx. 1-732-578-9777

www.mtc-usa.com Email: info@mtc-usa.com

For more information visit www.mtc-usa.com