## **Product Sheet**



## Catalog Number

## 49302-CX-75-010


**Union, Connector, Cross,** high pressure, 0.75 mm / 0.030 inch through hole with 10-32 screw threads. PEEK, round with an internal "X" flow pattern. ARE-Applied Research brand. 10/PK.



## Additional Info

This low dead volume, "Union" is used to connect PEEK or stainless steel tubing with the same connectors on both ends in the high pressures zones of HPLC systems. A "Cross (X)" internal flow pattern is often used to split flow into two or three directions or to combine two or three streams into one.

MicroSolv Technology Corporation 9158 Industrial Blvd. NE, Leland, NC 28451 tel. (732) 380-8900, fax (910) 769-9435 Email: customers@mtc-usa.com Website: www.mtc-usa.com

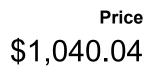

Printed On 2024/05/20