


Catalog Number

49681-TF-005

Inline Filter, for HPLC, female connections All-PEEK filter holder only. Uses disposable filters which are sold separately. ARE-Applied Research brand. 5/PK.

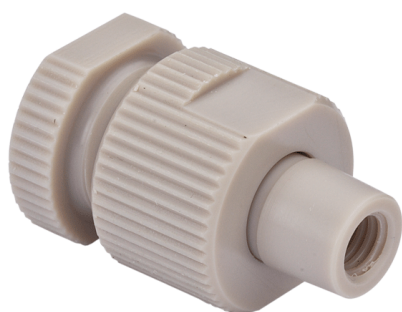

Pack of 5

Additional Info

Compatible with typical flow rate analytical columns up to 4,500 psi and 5 ml/minute. Also used as a union to connect tubing, but with a filter. *Not recommended for use with UHPLC pressure ranges.*

- Economical way to protect your columns and instruments.
- Re-usable filter holder.
- Easy to install after the pump, before the injector.
- Easy to replace the filters.

The universal connections (fits all standard HPLC threads and pilot lengths) can be installed with a wrench or hand tightened.

When these filters are in use, and pressure starts to rise, simply disconnect (after pressure in the column is at atmospheric), dispose of and replace with a suggested, new 0.188 inch, 0.5 um

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

tel. (732) 380-8900, fax (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com

Price

\$0.00

Printed On 2024/05/20