
## **Product Sheet**



## Catalog Number

## 9502C-S34-B

**Caps, Snap Top,** with 1.00mm thick, White Silicone Rubber / Red PTFE Pre-Slit (Y Shape) Septa, fitted in smooth, soft polyethylene, Blue Caps. For use with 11mm, 12x32mm Snap Top autosampler vials. MicroSolv

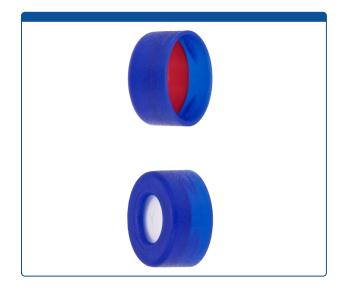

## Additional Info

Soft polyethylene caps are more pliable than other Snap Caps and are easier to place on the vial. Popular choice when using a large number of caps at one time. These caps may not seal as well as other Snap Caps and should be used with sample solvents with low volatility.

Y Shape Pre-Slit provides a greater target area for needles and are easier to penetrate and the 1.0mm Septa is thicker than most and offers improved sealing compared with other Pre-Slit Septa.

Pre-Slit Septa are normally recommended for Waters Alliance and Acquity instruments.

MicroSolv Technology Corporation 9158 Industrial Blvd. NE, Leland, NC 28451 tel. (732) 380-8900, fax (910) 769-9435 Email: customers@mtc-usa.com Website: www.mtc-usa.com



Printed On 2024/05/20