


# WINDOW MAKER™

Make detection windows for CE easily!



Reproducibly create a 2mm window in your CZE capillaries. Size of the window will be the same capillary to capillary and analyst to analyst. Different heating modules for alternate window sizes are sold separately.

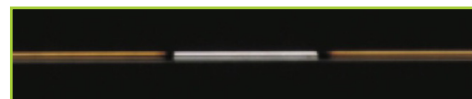

Example of window created by The Window Maker™

## Polyimide Coating Remover for HPCE

This unique unit is highly recommended for controlled stripping of polyimide from fused silica tubing creating a 2mm detection window that will not warp the glass during heating. With this controlled temperature the capillary remains intact as the burn area is supported on either side by polyimide. In just a few seconds you have a clean, window ready for detection with your CZE instruments. \*The Window Maker is not recommended for use on polyacrylamide, CEC capillaries or non CE based applications.

Cat. No.	Description	For more information visit <a href="http://www.mtc-usa.com">www.mtc-usa.com</a>
07200-S	Window Maker, MicroSolvCE™. Capillary coating stripper 110/120v includes a 2mm burn module.	
07201-S	Heating module for Window Maker. Replacement. Makes approximately a 2mm burn area.	
07204-S	Heating module accessory for Window Maker. Makes approximately a 4mm burn area.	
07207-S	Heating module accessory for Window Maker. Makes approximately a 7mm burn area.	

MANUFACTURED BY:

**MICROSOLV** TECHNOLOGY CORPORATION

9158 Industrial Blvd NE  
Leland, NC 28451

p: 1.732.380.8900  
f: 1.910.769.9435

customers@mtc-usa.com  
[www.Cogent-HPLC.com](http://www.Cogent-HPLC.com)