


Small Volume Inserts

Polypropylene for 2 mL, 12 x 32 mm Autosampler Vials

CONVENIENT • BIO-COMPATIBLE • LOW COST



These convenient limited volume inserts can be used with most autosampler vials that are 12 x 32 mm in dimension. Completely compatible with all MicroSolv 2 mL screw top vials or Snap Top vials these inserts are made from virgin, medical polypropylene and assure purity and inertness. Supplied with attached springs for centering and cushioning, these “all-purpose” inserts are an excellent choice when you are limited in the volume of sample you have and or want to maximize the sample injection without exposure to glass.

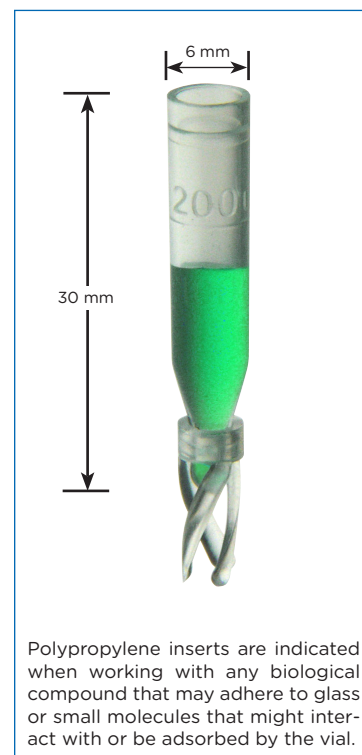

Polypropylene inserts are indicated when working with any biological compound that may adhere to glass or small molecules that might interact with or be adsorbed by the vial.

Cat. No.

9502S-02PP

Description

Small volume precision point polypropylene inserts 6 x 30 mm with attached plastic springs, with graduations 200 μ l, for 2 ml, 12 x 32 mm vials 100/pack. AQ™ Advanced Quality.

For more information visit www.mtc-usa.com

MANUFACTURED BY:

MICROSOLV TECHNOLOGY CORPORATION

9158 Industrial Blvd NE
Leland, NC 28451

p: 1.732.380.8900
f: 1.910.769.9435

customers@mtc-usa.com
www.Cogent-HPLC.com