



Are the Ends of the Pre Cut Stainless Steel Tubing Polished - FAQ

Date: 2-DECEMBER-2016 Last Updated: 3-NOVEMBER-2025


No, the ends of the pre-cut stainless-steel tubing are not polished. They are produced with a cutting device known as ECM (electro chemical machining) which produces a very flat, burr free edge which is better than polishing making them perfect for HPLC. See comparison of different types of cuts below.

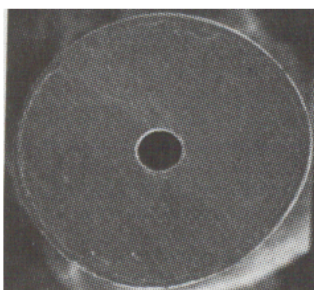

ECM-cut

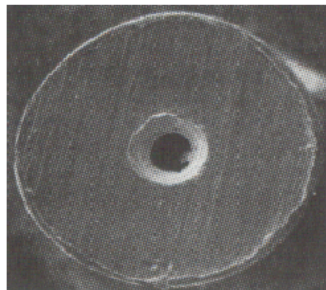

Mechanical cut

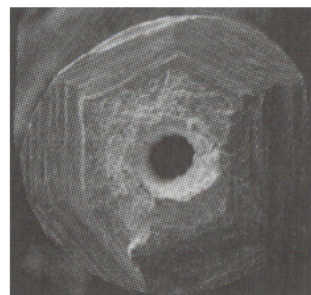

Manual cut

Click [HERE](#) for high pressure tubing ordering information and pictures.