
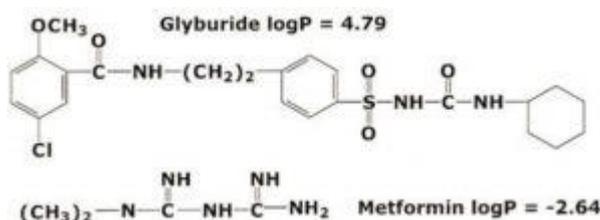
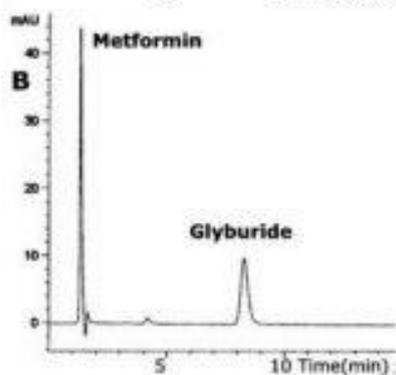


Glyburide and Metformin Analyzed with HPLC - AppNote

Separation of Highly Polar and Non-Polar Compounds in one Isocratic run

In this Method, the polar compound, Metformin, and the nonpolar compound Glyburide, can be retained on the same Stationary Phase (see A & B). Depending on the Mobile Phase composition either Metformin or Glyburide can be retained longer.

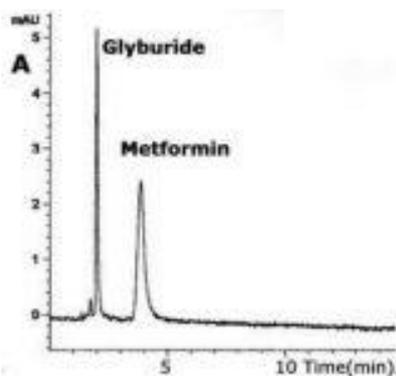

Method Conditions

Column: Cogent Bidentate C18™, 4 μm, 100 Å

Catalog No.: 40018-75P

Dimensions: 4.6 x 75 mm

Mobile Phase:

A: 15% DI Water / 85% Acetonitrile / 0.5% Formic Acid

B: 50% DI Water / 50% Acetonitrile / 0.5% Formic Acid

Injection vol.: 1 μ L
Flow rate: 0.5 mL / minute
Detection: UV @ 254 nm

Sample Preparation: Stock Solution: 100 μ g / μ L Glyburide and Metformin

Notes: Elution order was confirmed by LCMS, APCI+, with single ion monitoring Metformin (m/z 130) and Glyburide (m/z 369).



Attachment No 08 Glyburide and Metformin Analyzed with HPLC pdf 0.2 Mb [Download File](#)

Printed from the Chrom Resource Center
Copyright 2025, All Rights Apply
MicroSolv Technology Corporation
9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900
Fax: (910) 769-9435
Email: customers@mtc-usa.com
Website: www.mtc-usa.com