


Ticagrelor Analyzed with HPLC - AppNote

A Reproducible Method for Detection of this Blood Clot Preventative


Click [HERE](#) for Column Ordering Information.

A rapid, sensitive, and Reproducible Method has been developed for analysis of Ticagrelor. The data below, (*overlay of 5 chromatograms*) illustrates how the compound can be adequately Retained and detected using a simple Gradient in Reversed Phase HPLC. The Method demonstrates good Peak Shape and run-to-run Precision with RSD values less than 0.3%.

A Phenyl ring in the Column Stationary Phase provides beneficial π - π Interaction with the Analyte making possible the use of a very simple, Mass Spec friendly Mobile Phase with Formic Acid as an additive.

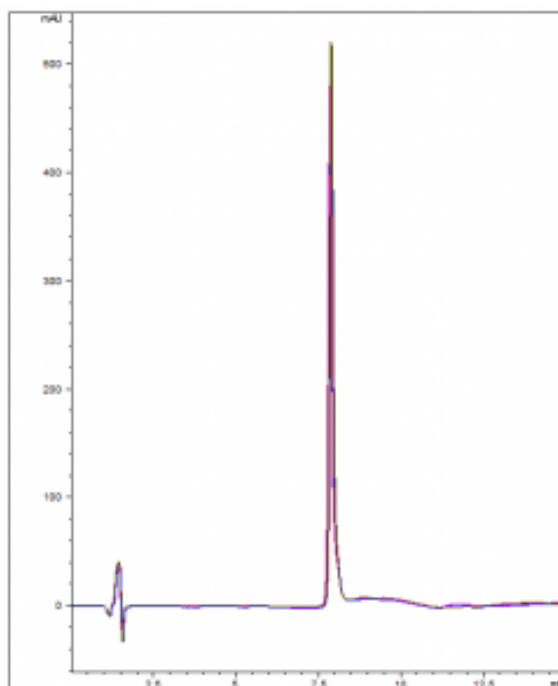

5 Injections of Ticagrelor Overlay

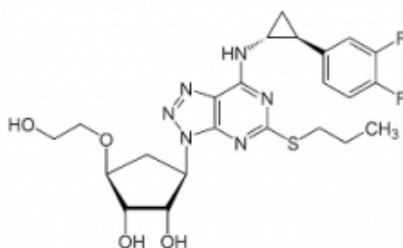

Ticagrelor

Method Conditions

Column: Cogent Phenyl Hydride™, 4µm, 100Å

Catalog No.: [69020-10P-2](#)

Dimensions: 2.1mm x 100mm

Mobile Phase:

A: DI Water with 0.1% Formic Acid (v/v)

B: Acetonitrile with 0.1% Formic Acid (v/v)

Gradient:

| Time (minutes) | %B |
|----------------|----|
| 0 | 40 |
| 1 | 40 |
| 5 | 85 |
| 6 | 85 |
| 7 | 40 |
| 10 | 40 |

Injection vol.: 1µL

Flow rate: 0.2mL / minute

Detection: UV @ 254nm

Sample Preparation: 2.0mg / mL Ticagrelor standard solution in Methanol.

Notes : Ticagrelor is an oral prescription indicated for treatment of acute ischemic stroke, coronary artery disease and acute coronary syndrome.



Printed from the Chrom Resource Center
Copyright 2025, All Rights Apply
MicroSolv Technology Corporation
9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900
Fax: (910) 769-9435
Email: customers@mtc-usa.com
Website: www.mtc-usa.com