


Sulfur Analyzed with HPLC - AppNote

Analyzing Elemental Sulfur Content

A simple and reproducible method has been developed for analysis of elemental Sulfur . The data below, (**an overlay of 5 chromatograms**) illustrates how the compound can be adequately retained in HPLC with good precision and peak shapes using this straightforward method.

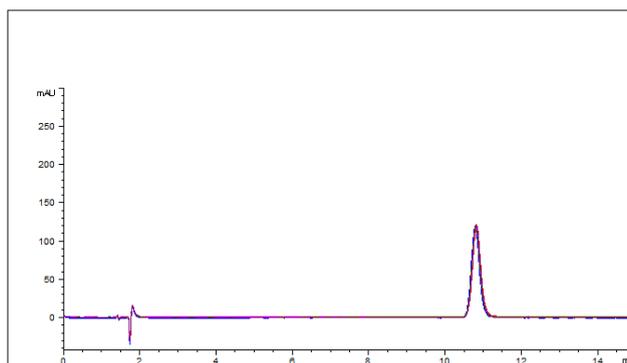

Peak: Sulfur

Method Conditions:

Column: Cogent Bidentate C18™, 4 µm, 100 Å.

Catalog No.: [40018-15P](#)

Dimensions: 4.6 x 150 mm

Mobile Phase: (90:10) Acetonitrile / DI water with 0.1% formic acid.

Injection vol.: 3 µL

Flow rate: 1.0 mL / minute

Detection: UV @ 263 nm

Sample Preparation: 0.5 mg / mL Sulfur (7704-34-9) in DCM

%RSD: <0.1%

t₀: 1.7 minutes

K': 5.2

Notes: This method could be used for determination of elemental Sulphur in soils, using DCM as the extracting solvent.

Note 2: Capacity was determined using the following equation: $k = (t_R - t_0) / t_0$

t_R = Retention time of an analyte peak

t_0 = Retention time of non-retained peak



Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com