
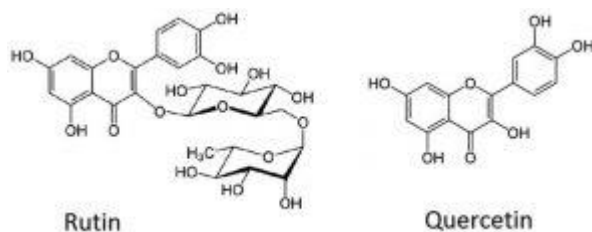


## Quercetin Rutin in Brazilian Tree Bark Analyzed with LCMS - AppNote

### Identification of Quercetin, Rutin, and "Compound 1" in a Brazilian Tree Bark

Over one quarter of natural medicines have been discovered in rainforests, and therefore analytical study of the indigenous flora and fauna is necessary. The Amazonian rainforest in Brazil is the most biodiverse region in the world.

A Brazilian tree bark extract (*Brownea Grandiceps*) is analyzed by LCMS and three beneficial compounds were identified in the Extracted Ion Chromatograms (EIC)s.

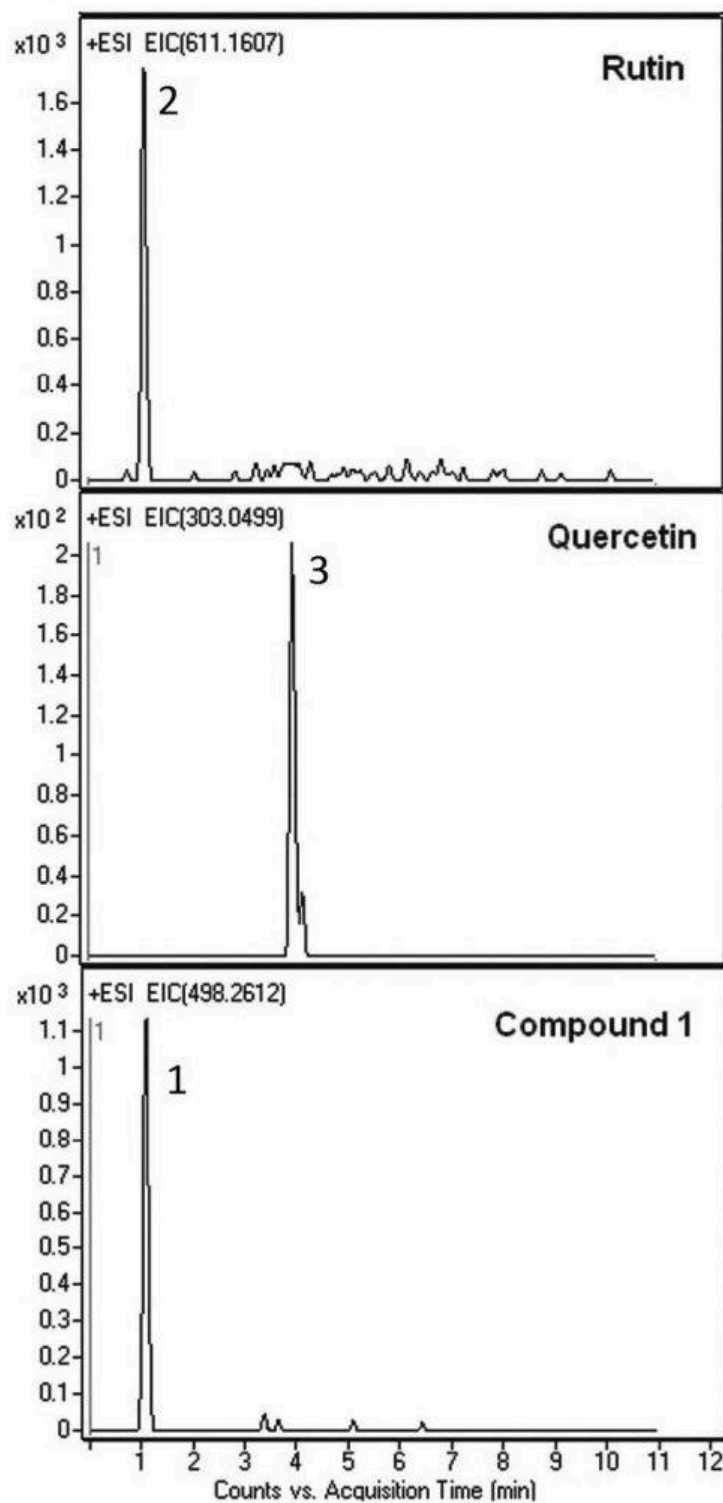

### Peaks:

1. (6-Beta-O-2',3'-Dihydrocinamonyl-12-Hydroxy-(13) 15-en-16,12-Olide-18-Cassaneoic Acid)  $m/z$  498.2612  $[M + H]^+$
2. Rutin  $m/z$  611.1607  $[M + H]^+$
3. Quercetin  $m/z$  303.0499  $[M + H]^+$

## Method Conditions

**Column:** Cogent Diamond Hydride™, 4µm, 100Å

**Catalog No.:** [70000-15P-2](#)

**Dimensions:** 2.1 x 150mm

**Mobile Phase:**

A: DI Water / 0.1% Formic Acid

B: Acetonitrile / 0.1% Formic Acid (v/v)

**Gradient:**

Time (minutes)	%B
0	80
1	10
5	10
6	80

**Post Time:** 3 minutes

**Injection vol.:** 1µL

**Flow rate:** 0.4mL / minute

**Detection:** ESI – POS - Agilent 6210 MSD TOF Mass Spectrometer

**Sample Preparation:** Two pieces of bark were boiled for 5 minutes in DI Water and then filtered with a 0.45µm Nylon Syringe Filter (MICROSOLV Tech Corp.).

**t<sub>o</sub>:** 0.9 minutes

---

*Note: Preliminary research suggests that Quercetin may have Antiviral, Anti-Cancer, and Anti-Inflammatory properties. Likewise, Rutin has been reputed to have health benefits as well such as Anti-Oxidant properties, lowered risk of heart attack or stroke, and others.*

---



**Attachment No 320 Brazilian Tree Bark Analyzed with LCMS pdf 0.3 Mb** [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)