
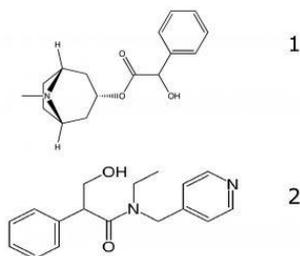


Tropicamide and Homatropine Analyzed by HPLC - AppNote

Separation of Anticholinergic Mydriatic Agents

In this Method, Tropicamide and Homatropine are separated with a simple Gradient. Homatropine is significantly hydrophobic yet still retains well.

Furthermore, the Separation between these two solutes is good, illustrating the Selectivity capabilities of this Method. Data is shown for two lots of Columns in order to demonstrate the robustness and precision of this Method.

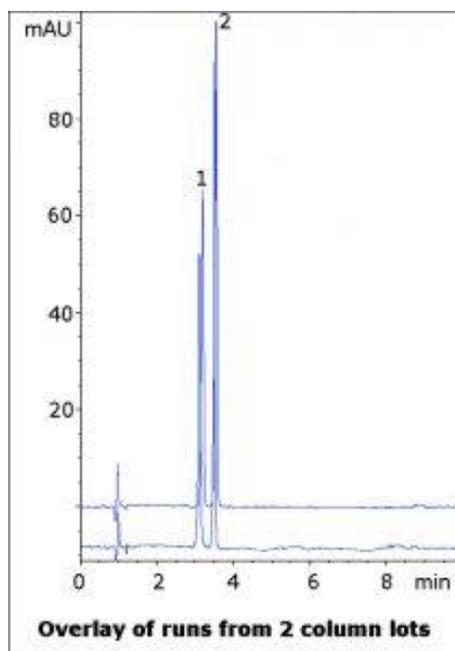

Peaks: 1. Homatropine, 2. Tropicamide

Method Conditions

Column: Cogent Diamond Hydride™, 4 μm, 100 Å

Catalog No.: 70000-7.5P

Dimensions: 4.6 x 75 mm

Mobile Phase:

- A: DI Water / 0.1% Trifluoroacetic Acid (v/v)
- B: Acetonitrile / 0.1% Trifluoroacetic Acid (v/v)

Gradient:

Time (minutes)	%B
0	95
1	95
6	50
7	95

Injection vol.: 2 µL**Flow rate:** 1.0 mL / minute**Detection:** UV @ 220 nm

Sample Preparation: 0.1 mg / mL Homatropine and 0.01 mg / mL Tropicamide reference standards in diluent of 50:50 Solvent A / Solvent B. Peak identities were confirmed with individual standards.

t_o: 0.9 minutes

Note: Homatropine and Tropicamide are mydriatic (pupil-dilating) agents used in ophthalmoscopic examinations. Trade names include Homatropaire® for Homatropine and Mydral® for Tropicamide.



Attachment No 219 Tropicamide and Homatropine Analyzed by HPLC pdf 0.4 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.comWebsite: www.mtc-usa.com