
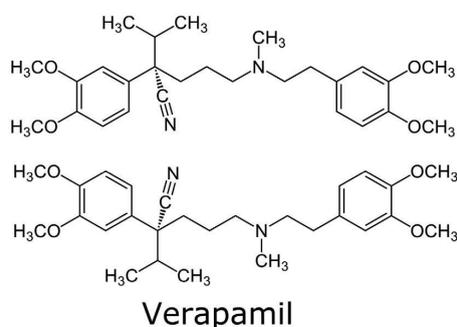


Verapamil and Impurity Analyzed with HPLC - AppNote

Alternative to USP Method with Improved Peak Efficiency

The USP impurities Method for Verapamil uses an Acetate Buffer with the ion pair agent 2-Aminoheptane to reduce Peak Tailing of the API and its Impurity. This Application Note uses an LCMS compatible Mobile Phase and produces Symmetrical Peaks for both Compounds. Excellent Resolution was also obtained between the two Peaks, far exceeding the USP System Suitability of $R_s \geq 1.5$.

This Method is a significant improvement over the USP monograph and illustrates the considerable potential for analyses of Amine containing solutes.

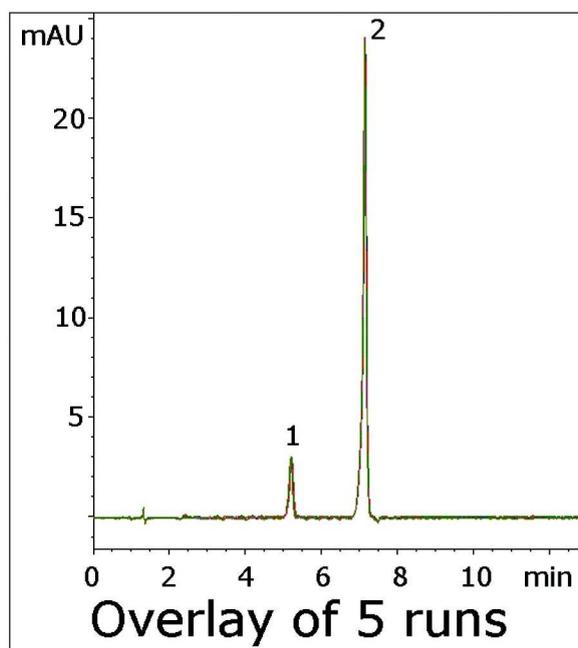

Peaks: 1. Verapamil Related Compound B, 2. Verapamil

Method Conditions

Column: Cogent Diamond Hydride™, 4 μm, 100 Å

Catalog No.: 70000-7.5P

Dimensions: 4.6 x 75 mm

Mobile Phase:

- A: DI Water / 10mM Ammonium Formate
- B: 95:5 Acetonitrile / 10mM Ammonium Formate (v/v)

Gradient:

Time (minutes)	%B
0	100
9	85
10	100

Post Time: 4 minutes

Injection vol.: 1 µL

Flow rate: 1.0 mL / minute

Detection: UV @ 278 nm

Sample Preparation: Sample mixture contained 0.1 mg / mL Verapamil and 0.02 mg / mL Verapamil Related Compound B in a diluent of 50:50 Solvent A / Solvent B. Peak identities were confirmed with individual Standards.

t_o : 0.9 minutes

Note: Verapamil is an L-type calcium channel blocker used in the treatment of Hypertension, Angina Pectoris, and Cardiac Arrhythmia.



Attachment No 180 Verapamil and Impurity Analyzed with HPLC pdf 0.5 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com