
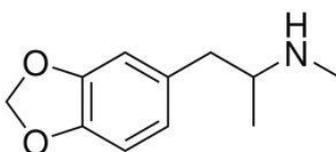


MDMA Analysis in Plasma Samples with LCMS - AppNote

Methylenedioxyamphetamine Analyzed with LC-MS

Click [HERE](#) for Column Ordering Information.

Under the described conditions, MDMA was retained and eluted as a Symmetrical Peak. The Sensitivity of the Method is very good and comparable to that reported with GCMS Detection [1]. Matrix effects were of minor extent and reproducible and hence should not compromise Quantification. The Method can be used for Forensic Research and Clinical Analysis.

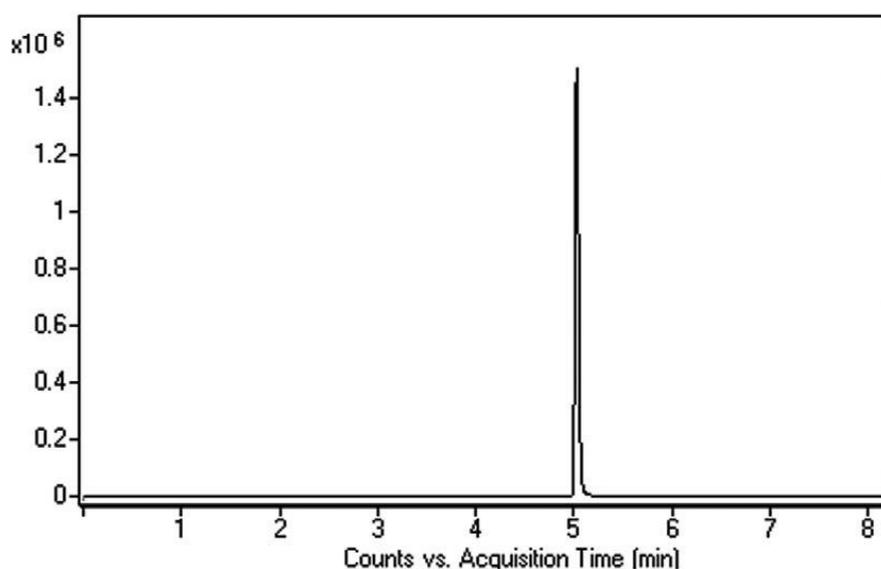

Peak:

(±)-3,4-Methylenedioxyamphetamine, m/z 194.1176 [M+H]⁺

Method Conditions

Column: Cogent Phenyl Hydride™, 4 μm, 100 Å

Catalog No.: [69020-05P-2](#)

Dimensions: 2.1 x 50 mm

Mobile Phase:

A: DI Water / 0.1% Formic Acid (v/v)

B: Acetonitrile / 0.1% Formic Acid (v/v)

Gradient:

Time (minutes)	%B
0	10
1	10
6	90
7	10

Post Time: 3 minutes

Flow rate: 0.4 mL / minute

Injection vol.: 1 µL

Sample Preparation: 50 µl of Acetonitrile was mixed with 50µl of plasma for protein precipitation. The samples were centrifuged (16000×g for 15 minutes) , and the supernatant was filtered through a 0.45 µm Nylon Syringe Filter (MICROSOLV Tech Corp.) and transferred to autosampler vials for injection.

Detection: ESI – POS - Agilent 6210 MSD TOF Mass Spectrometer

t_o : 0.9 minutes

Note: The Amphetamine derivative 3,4-methylenedioxyamphetamine (MDMA), known also as Molly or Ecstasy, is often used or abused as a recreational drug. Because of a reported high inter-individual difference of its toxicity, sensitive analytical methods are needed. A urine test is a standard method to investigate drug abuse but the method has a very low diagnostic sensitivity and makes testing in plasma much more suitable.

Reference: [1]. R. Kikura, Y. Nakahara, T. Mieczkowski, F. Tagliaro, *Forensic Sci. Int.* 84 (1997) 165–177.



Attachment No 263 Analysis of MDMA in Plasma Samples with LCMS pdf 0.2 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com