

## Formoterol Fumarate Analyzed with HPLC - AppNote


### Forced Degradation Methods


Click [HERE](#) for Column Ordering Information.

Forced degradation studies are useful for developing Stability Indicating Methods of Pharmaceuticals. In this Application Note Separation is obtained between the API and a Degradant formed under two different, strong acidic conditions. The non-degraded Sample produced a sharp, Symmetrical Peak for the API shown in *Figure A*.

In *Figure B*, strong acid conditions without heating produced a degradant while the API Peak was still present. In *Figure C*, acid degradation with heating completely converted the API into the degradant, as the API Peak was no longer observed.

Five Chromatograms are overlaid and presented below to show Method Robustness and Precision.

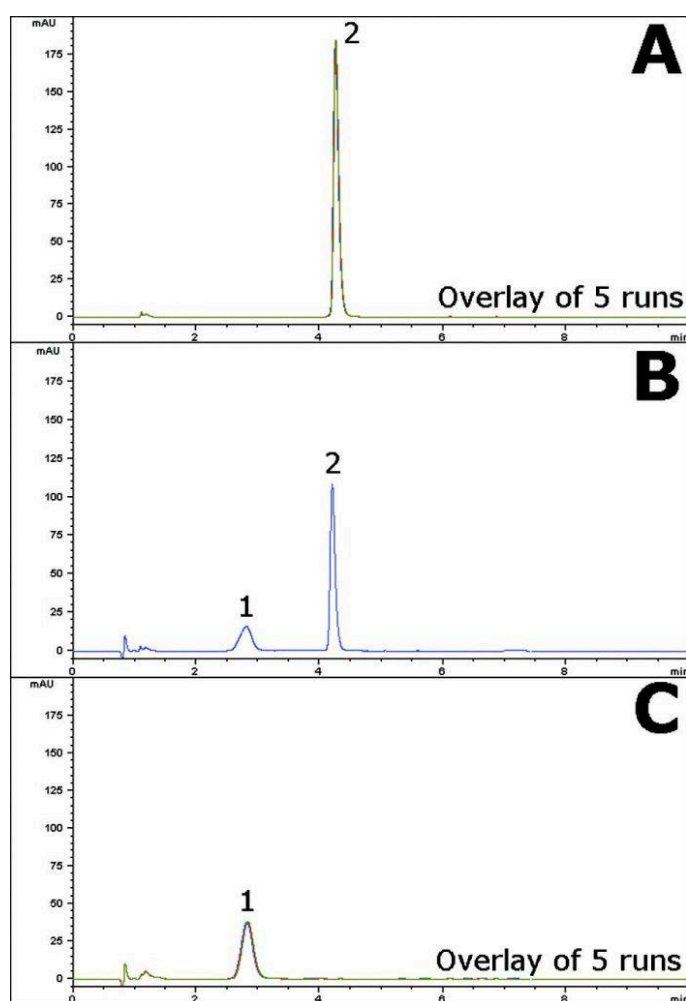

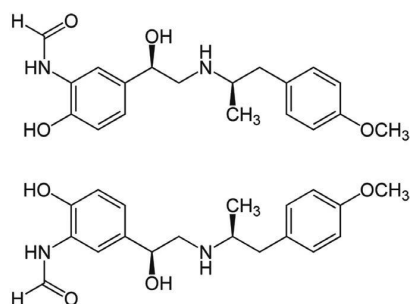

### Peaks:

1. Degradant
2. Formoterol Showing Enantiomers

### Method Conditions

**Column:** Cogent Phenyl Hydride™, 4μm, 100Å

**Catalog No.:** 69020-7.5P

**Dimensions:** 4.6 x 75mm

**Mobile Phase:**

A: DI Water / 0.1% Trifluoroacetic Acid (TFA) v/v

B: Acetonitrile / 0.1% Trifluoroacetic Acid (TFA) v/v

**Gradient:**

Time (minutes)	%B
0	10
1	10
7	40
8	10

**Post Time:** 2 minutes

**Injection vol.:** 10μL

**Flow rate:** 1.0mL / minute

**Detection:** UV @ 282nm

**Sample Preparation:**

*Figure A:* Non-degraded Sample: 0.1mg / mL Formoterol Fumarate in 50:50 Solvent A / Solvent B Diluent.

*Figure B:* Acid degradation of the Sample, with no heating: 0.1mg / mL Formoterol Fumarate in 50:50 1N HCL / Acetonitrile Diluent.

*Figure C:* Acid degradation of the Sample, with heating: 0.1mg / mL Formoterol Fumarate in 50:50 1N HCL / Acetonitrile Diluent after 30 minutes of heating at 85°C.

**t<sub>0</sub> :** 0.9 minutes

**Note:** Formoterol is used a long-acting β<sub>2</sub>-antagonist used in treatment of conditions such as asthma and obstructive pulmonary disease. It is marketed under several trade names such as Foradil® and Oxeze®.

**Attachment No 178 Formoterol Fumarate pdf** 0.4 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

**MicroSolv Technology Corporation**

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: [customers@mtc-usa.com](mailto:customers@mtc-usa.com)

Website: [www.mtc-usa.com](http://www.mtc-usa.com)