
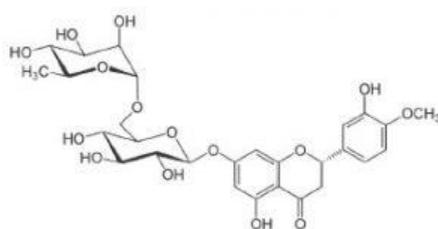


Hesperidin Analyzed with HPLC - AppNote

LC-MS compatible method

This method for analysis of Hesperidin produces an excellent analyte peak and shows separation from what appears to be two impurity peaks. The gradient method is LC-MS compatible and can be applied to many types of Hesperidin-containing samples such as citrus fruit extracts. The overlay of five runs shown in the figure illustrates the run-to-run reproducibility of the data.

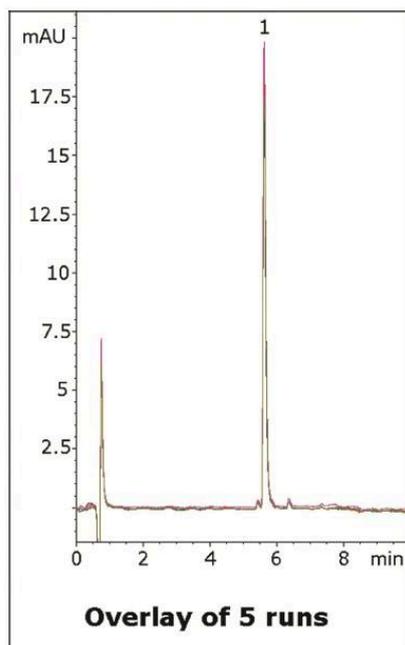

Peak:
Hesperidin

Method Conditions

Column: Cogent Bidentate C18 2.0, 2.2 μm , 120 \AA

Catalog No.: [40218-05P-2](#)

Dimensions: 2.1 x 50 mm

Mobile Phase:

A: DI Water / 0.1% Formic Acid (v/v)

B: Acetonitrile / 0.1% Formic Acid (v/v)

Gradient:

Time (minutes)	%B
0	10
1	10
9	70
10	10

Post Time: 5 minutes**Injection vol.:** 1 µL**Flow rate:** 0.3 mL/minute**Detection:** UV @ 285 nm**Sample:** 100 ppm Hesperidin reference standard in 1:1 DMSO: Methanol diluent**t₀:** 0.6 minutes

Note: Hesperidin is a naturally occurring glycoside polyphenol that is thought to have antioxidant properties. Studies suggest it may have a number of pharmaceutical applications due to possible anti-inflammatory, anti-cancer, and cholesterol and blood pressure lowering effects. Its name comes from the word "hesperidium", which is the type of fruit produced by citrus trees.

**Attachment No 289 Hesperidin Analyzed by HPLC pdf** 0.3 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.comWebsite: www.mtc-usa.com