


HPLC Method Transferred to Near UHPLC - AppNote

Date: 10-JUNE-2014 Last Updated: 7-FEBRUARY-2026

Use of near UHPLC column for improved results

This challenging compound has several amine groups and can be problematic to analyze by HPLC. Here a good peak shape is obtained using the Cogent Diamond Hydride 2.0 Column. The efficiency is higher on the 2.2 μ m phase compared to a standard 4 μ m Column, leading to increased sensitivity. The method conditions are compatible with LCMS as well.

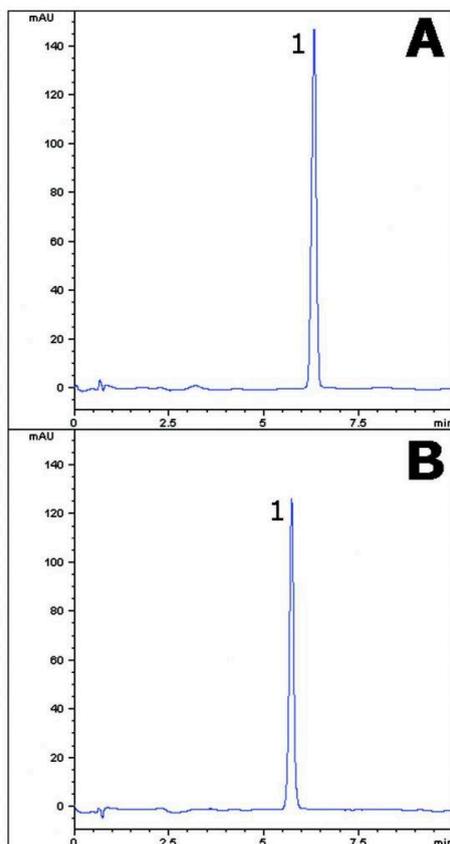

Tizanidine

Peak: Tizanidine HCl

Method Conditions

Columns:

- Fig. A: Cogent Diamond Hydride 2.0™, 120Å
- Fig. B: Cogent Diamond Hydride™, 4µm, 100Å

Catalog Nos.:

- Fig. A: [70200-05P-2](#);
- Fig. B: [70000-05P-2](#)

Dimensions:

- Fig. A: 2.1 x 50 mm
- Fig. B: 2.1 x 50 mm

Mobile Phase:

- A: DI Water / 0.1% Formic Acid (v/v)
- B: Acetonitrile / 0.1% Formic Acid (v/v)

Gradient:

Time (minutes)	%B
0	95
1	95
6	40
7	95

Post time: 3 minutes

Injection vol.: 0.2 µL

Flow rate: 0.3 mL/minutes

Detection: UV @ 230 nm

Sample: 4 mg strength Tizanidine HCl tablet was ground and weighed in a 10 mL Volumetric Flask. A portion of 50/50 Solvent A / Solvent B diluent was added and the flask was sonicated 10 minutes. It was then diluted to mark and filtered with a 0.45 µm Nylon Syringe Filter (MICROSOLV Tech Corp.).

t_o: 0.6 minutes

Note: Tizanidine is a centrally acting A2-Adren- Ergic agonist used to treat spasms, cramping, tightness of muscles, and related conditions. It is available under the trade name Zanaflex® as well as generic versions.

Attachment Tizanidine HCl Method Transfer pdf [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com