



Vancomycin HPLC Method Transferred to Near UHPLC - AppNote

Benefits of using 2.2 μm column for Vancomycin analysis

The highly polar antibiotic Vancomycin can be readily analyzed with either the standard 4 μm Cogent Diamond Hydride or the Cogent Diamond Hydride 2.0 column. Notably higher efficiency is obtained on the smaller particle size column.

This Method is easy to perform and is LCMS compatible. Three runs were performed on each column in order to demonstrate consistency and show the benefits of using the smaller particle column.

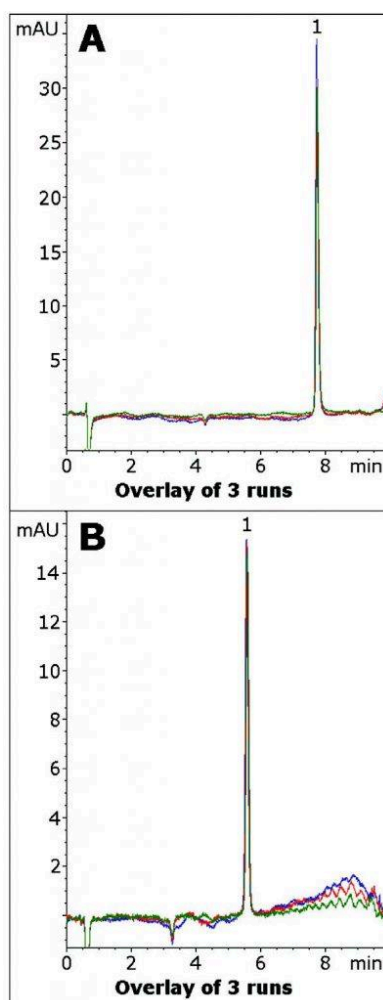

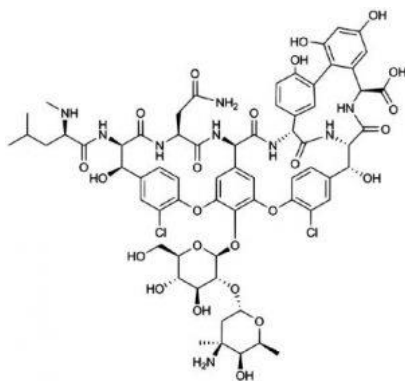

Peak:
Vancomycin

Method Conditions

Columns:

Fig. A: Cogent Diamond Hydride 2.0™, 2.2µm, 120Å

Fig. B: Cogent Diamond Hydride™, 4µm, 100Å

Catalog Nos.:

Fig. A: [70200-05P-2](#);

Fig. B: [70000-05P-3](#)

Dimensions:

Fig. A: 2.1 x 50 mm

Fig. B: 3.0 x 50 mm

Mobile Phase:

A: DI Water / 0.1% Formic Acid (v/v)

B: Acetonitrile / 0.1% Formic Acid (v/v)

Time (minutes)	%B
0	70
0.4	70
7.4	10
8.4	70

Post time: 5 minutes

Injection vol.: 0.2 µL

Flow rate:

Fig. A: 0.29 mL / minutes

Fig. B: 0.50 mL / minutes

Detection: UV @ 210 nm

Sample:

Stock Solution: 1 mg / mL Vancomycin HCl in 50 / 50 solvent A / solvent B diluent. The solution was filtered through a 0.45µm Nylon Syringe Filter (MICROSOLV Tech Corp.).

Working Solution: Stock solution was diluted 1:100 with 50/50 solvent A/ solvent B mixture.

t₀: 0.4 minutes

Note: Vancomycin is a glycosylated non-ribosomal Peptide antibiotic used to treat Colitis. Vancomycin is often used as a drug of last resort when other antibiotics are rendered ineffective due to developed resistance of bacteria. It is a natural product isolated from Amycolatopsis Orientalis.



Attachment Vancomycin Method Transfer pdf [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com