


Cogent Diamond Hydride 2.0 HPLC column retention time precision - AppNote

This graph is from a study investigating the day to day retention time precision with two test solutes on the Cogent Diamond Hydride 2.0™ HPLC column

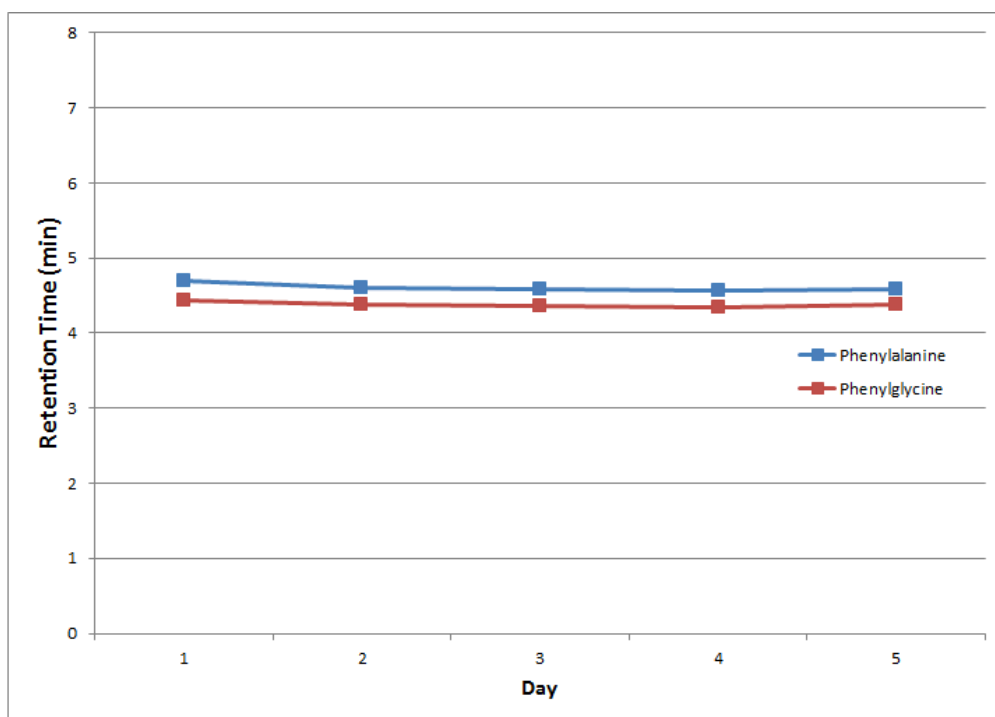

Peaks:

Blue Trace: Phenylalanine

Red Trace: Phenylglycine

Many factors can affect column robustness and lifetime, so it is not possible to define specific values for those parameters.

Mobile Phase:

A: DI Water / 0.1% Formic Acid

B: Acetonitrile / 0.1% Formic Acid

Flow Rate: 0.3mL / minute

Injection Volume: 0.2uL

Detection: UV 254nm

Gradient:

Time (minutes)	%B
0	90

6	20
7	20
8	90



[Cogent Diamond Hydride HPLC Column Ordering Information](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com