
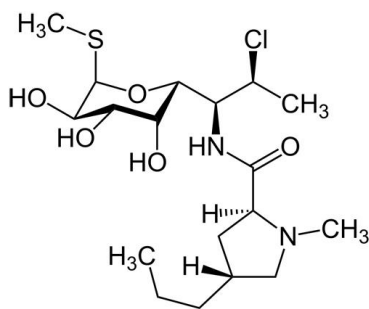


Clindamycin Antibiotic Tablet Analyzed with HPLC - AppNote

Separation of API from Excipients

The Stability and Robustness of this Method makes it an excellent choice for QC of this semi-synthetic antibiotic which is structurally related to Lincomycin. The Chromatogram shown below of Clindamycin is from a complex Pharmaceutical Formulation but note the excellent Peak Shape for it and the Separation from the Excipients.

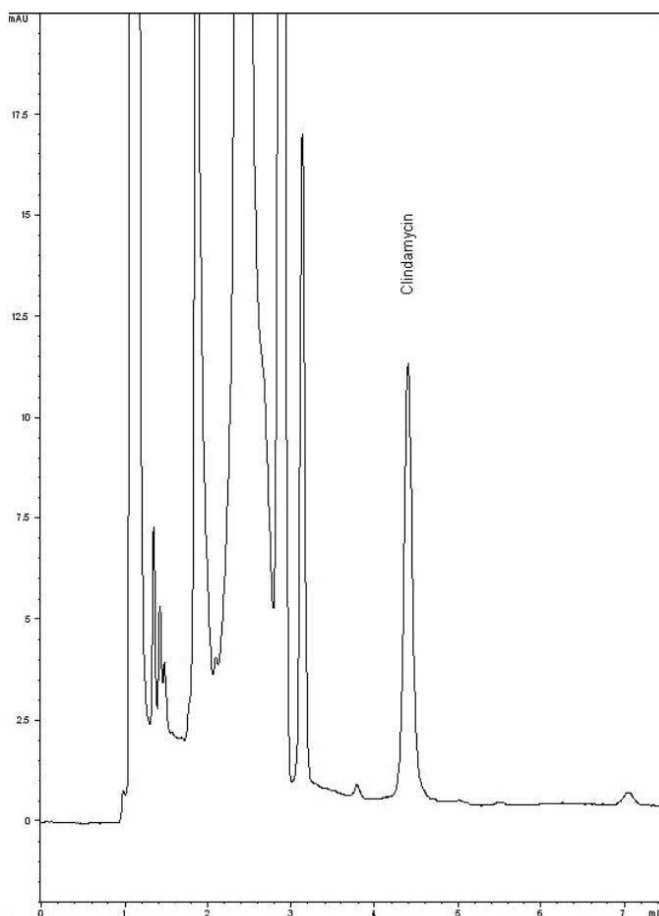

Peak:
Clindamycin

Method Conditions

Column: Cogent Bidentate C18™, 4μm, 100Å

Catalog No.: [40008-15P](#)

Dimensions: 4.6 x 150mm

Mobile Phase: 40:60 Acetonitrile / 10mM Phosphate Buffer pH 7.5

Temperature: 25°C

Injection vol.: 20μL

Flow rate: 1.0mL / minute

Detection: UV @ 338nm

Notes: Normally used to treat infections with anaerobic bacteria, this drug is indicated for use against Protozoal diseases such as Malaria. Also used for Acne can be useful against some MRSA infections. Clindamycin is marketed under various trade names including Dalacin, Cleocin, Evoclin and Duac as well in generic form. Molecular Weight of the Clindamycin is 424.99 with a pKa of 7.6 and is soluble in Water, Pyridine, Ethanol and DMF.

If you use the Phosphate Buffer on the Bidentate C8 Column, you should dedicate the Column to use with Phosphates in the Mobile Phase.



Attachment No 71Clindamycin Antibiotic Tablet Analyzed with HPLC pdf 0.2 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.com

Website: www.mtc-usa.com