


Gly-Tyr, Val-Tyr-Val, Methionine Enkephalin, Angiotensin II, Leucine Enkephalin Peptide Mixture Separated by HPLC - AppNote

Precise and Fast Equilibration Time for this Peptide Sample

HPLC is a technique used for the characterization of Peptides and Reversed Phase HPLC is employed for the initial analysis and final, large scale purification. The first step of the production of Synthetic Peptides usually involves an initial Separation of these compounds from its mixture, on an analytical scale, then purification and collection of the target Peptide follows.

The Figure shown here presents the use of Reversed Phase HPLC Method in the Separation of a five Peptide Mixture. The 300Å pore size of the sorbent is ideal for separation of these Peptides chosen.

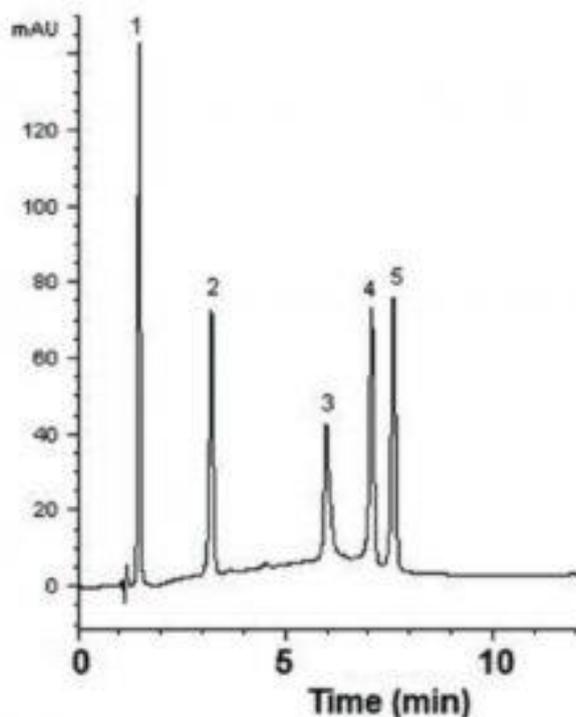

Peaks:

1. Gly-Tyr, 2. Val-Tyr-Val, 3. Methionine Enkephalin, 4. Angiotensin II, 5. Leucine Enkephalin

Method Conditions

Column: Cogent Bidentate C8 300™, 5 µm, 300 Å

Catalog No.: 40008-75P-3M

Dimensions: 4.6 x 75 mm

Mobile Phase:

A: DI Water / 0.1% Trifluoroacetic Acid (TFA)

B: Acetonitrile / 0.1% Trifluoroacetic Acid (TFA)

Gradient:

Time (minutes)	%B
0	9
5	21
20	27
21	9

Post Time: 5 minutes**Flow rate:** 1.0 mL / minute**Detection:** UV @ 214 nm

Notes: Peptides are distinguished from proteins on the basis of the number of amino acid residues. Generally, this number is about 50 residues or fewer. Because of their smaller size, Peptides do not have the same degree of complexity that proteins do.



Attachment No 91 Peptide Mixture Separation by HPLC pdf 0.2 Mb [Download File](#)

Printed from the Chrom Resource Center

Copyright 2025, All Rights Apply

MicroSolv Technology Corporation

9158 Industrial Blvd. NE, Leland, NC 28451

Tel: (732) 380-8900

Fax: (910) 769-9435

Email: customers@mtc-usa.comWebsite: www.mtc-usa.com